



UVB Narrowband PL-L/PL-S – most effective phototherapy lamp plus design freedom

UV-B Narrowband PL-L/PL-S

More than 400 independent clinical studies have proven that the UVB Narrowband PL-L/PL-S lamps are safer and more effective than any other lamps in their class. The PL-L/PL-S versions provide additional design freedom for system manufacturers since these are space-saving compact, single-ended lamps. Further flexibility is assured since they use the same lamp caps as general lighting lamps, they also use the same universal ballasts. Their best-in-class treatment effectiveness is a result of their very narrow waveband emission: between 305 and 315 nm with a peak at 311 nm. This is the most efficacious waveband for the treatment of psoriasis, and is therefore much more focused. As a result, exposure times are much shorter and this in turn leads to a reduction of side effects such as reddening of the skin and itching. All of this makes them ideal for phototherapy treatment of diseases such as psoriasis and vitiligo. What's more, because the overall dosage of this narrowband radiation can be closely controlled, these lamps are suitable for home therapy.

Avantages

- Minimum side effects like redness, itching and burns
- Shorter period of exposure and less erythematous radiation than conventional UVB lamps
- Optimal therapeutic effect with minimum side effects
- Proven to be most effective on the skin

UV-B Narrowband PL-L/PL-S

Fonctions

- Emission peak at 311 nm
- Narrowband
- Special developed phosphor and glass
- World wide tested in more than 400 clinical tests

Application

- Psoriasis, Vitiligo

Mises en garde et sécurité


- Un bris de lampe est très peu susceptible d'influer sur votre santé. Si une lampe se brise, aérez la pièce pendant 30 minutes et retirez les morceaux, de préférence avec des gants. Placez-les dans un sac en plastique scellé et portez-le à votre site de gestion des déchets en vue de son recyclage. Évitez d'utiliser un aspirateur.

Versions



XPPR XUVBPLS 0001

Schéma dimensionnel

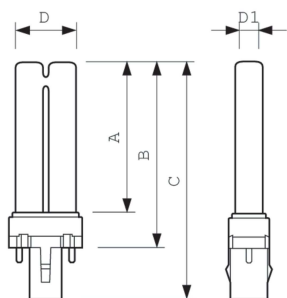

Product	D1 (max)	D (max)	A (max)	B (max)	C (max)
PL-S 9W/01/2P 1CT/6X10BOX	13 mm	28 mm	129 mm	144,5 mm	167 mm

UV-B Narrowband PL-L/PL-S

